



ELECTRODES & ACCESSORIES

●pH ●mV(ORP) ●ION ●Conductivity ●Dissolved Oxygen



CONTENTS

pH METER and ELECTRODE COMBINATION TABLE	1
pH ELECTRODES (3-in-1 ELECTRODES)	2
pH ELECTRODES (ISFET (Semiconductor electrode))	4
pH ELECTRODES (GLASS ELECTRODES(G), REFERENCE ELECTRODES(R))	5
pH ELECTRODES (COMBINATION), Temperature Compensation Electrode, METALLIC ELECTRODES (FOR ORP MEASUREMENT)	6
ION ELECTRODES	7
CONDUCTIVITY ELECTRODE CELLS	10
DISSOLVED OXYGEN(DO) ELECTRODE & TIPS	11
ACCESSORIES	12
STANDARD SOLUTIONS, INTERNAL SOLUTION for REFERENCE ELECTRODE & CLEANING SOLUTIONS	14
ACCESSORIES for U-50, U-20XD, U-10, W-20XD SERIES & INTERNAL SOLUTION for REFERENCE ELECTRODE	15
ELECTRODES & ACCESSORIES for TWIN/CARDY	16
pH Electrode Selection Guide	18
HORIBA WATER QUALITY ANALYZER LINEUP	20

pH ELECTRODES **METALLIC ELECTRODES** **ION ELECTRODES** **CONDUCTIVITY ELECTRODE CELLS** **DO ELECTRODES** **ACCESSORIES**

■ Applicable Product Models

Benchtop pH/Water Quality Analyzer	F-70 Series, F-50 Series, F-20/F-20 II Series, F-10 Series, M-10 Series
Portable pH Meter	F-70 Series, D-50 Series, D-20 Series, D-10 Series
Compact pH Meter	B-211/212/213/711/712/713
Benchtop Conductivity Meter	DS-70 Series, DS-50 Series, DS-10 Series
Portable Conductivity Meter	ES-71, ES-51, ES-10 Series
Compact Conductivity Meter	B-173/771
Compact Ion Meter	B-341/342/343/721/722/731/ 741/742/743/751, C-121/122/131/141
Portable DO Meter	OM-71, OM-51, OM-10 Series
Portable Water Quality Monitoring System	U-50, U-20XD, U-10 Series

pH METER and ELECTRODE COMBINATION TABLE

Type	pH					ORP	ION			Conductivity Electrode Cells	Dissolved Oxygen Electrode
	3-in-1 Electrode	Combination Electrode	ISFET Electrode	Single Electrode*1	Reference Electrode	3-in-1 Combination Electrode	Combination Electrode	Single Electrode*1			
	9615-10D	6069-10C	0030-10D	1066A-10C	2060A-10T	9300-10D	6560-10C	8001-10C	8011-10C	9382-10D	9520-10D
	9625-10D	6261-10C	0040-10D	1076A-10C	2565A-10T		6561-10C	8002-10C	1512A-10C	3551-10D	9551-20D
	9618-10D						5002A-10C	8003-10C	8201-10C	3552-10D	9551-100D
	9681-10D						6581-10C	8004-10C	8202-10C	3553-10D	
	9680-10D						6582-10C	8005-10C	8203-10C	3561-10D	
	6367-10D						6583-10C	8006-10C		3562-10D	
	6377-10D							8007-10C		3573-10C	
	6252-10D							8008-10C		3574-10C	
	9600-10D							8009-10C			
	9631-10D							8010-10C			
	9632-10D										
F-71, F-51・52	○	○	○	○	○	○	×	×		×	×
F-72・73, F-53	○	○	○	○	○	○	○	○		×	×
F-54	○	○	○	○	○	○	×	×		○	×
F-74・74BW, F-55	○	○	○	○	○	○	○	○		○	×
D-71, D-51, D-21	○	○	○	×	×	×	×	×		×	×
D-72, D-52, D-22	○	○	○	×	×	○	×	×		×	×
D-73, D-53, D-23	○	○	○	×	×	○	○	×		×	×
D-74, D-54, D-24	○	○	○	×	×	○	×	×		○	×
D-75, D-55, D-25	○	○	○	×	×	○	×	×		×	○
F-21・22・21Ⅱ・22Ⅱ	○	○	○	○	○	○	×	×		×	×
F-22C・22ⅡC	○	○	○	○	○	○	×	×		×	×
F-23・24・23Ⅱ・24Ⅱ	○	○	○	○	○	○	○	○		×	×
F-23C・24C・23ⅡC・24ⅡC	○	○	○	○	○	○	○	○		×	×
M-11, F-11・12	○	○*2	○	○*2	○*2	×	×	×		×	×
M-12・13, F-13・14・15・16	○	○	○	○	○	○	×	×		×	×
D-11・12	○	○*2	○	×	×	×	×	×		×	×
D-13・14	○	○*2	○	×	×	○	×	×		×	×

○: Not applicable ×: Applicable *1: Reference electrode required for measurement *2: Temperature compensation electrode (4163-10T) required for measurement

Electrode connector and lead wire length:

10 of -10C, -10T, or -10D in the last part of each type shows that the lead wire length is 1.0m. C, T, and D denote connector types for the main unit. The connector type suited for the main unit should be selected.

Only D type connector can be used for the D-20, D-50, D-70 series. C, T, D type connectors can be used for all of the F series and M-series.

<Reference>

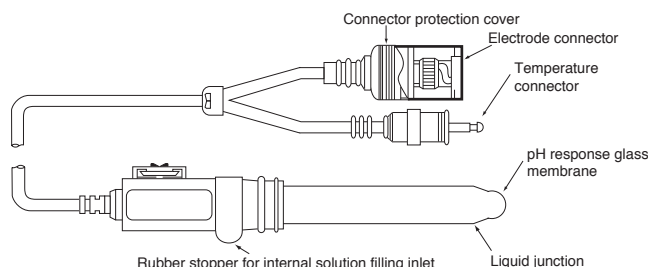
The liquid junction is the section where the liquid inside the reference electrode comes in contact with the sample liquid. Several junction types are available (ceramic, sleeve, etc.), to meet the requirements of specific samples or applications.

Liquid junction type	Features
Ceramic	A broad range of pH measurements. (Please note that samples of high viscosity may cause clogging.)
Movable sleeve	The larger liquid junction area is ideal for samples of high liquid junction potential, such as those with (1) high viscosity, (2) high salt concentration, or (3) low ionic strength. The liquid junction is easy to clean. High internal solution outflow volume.
Fixed sleeve	The large liquid junction area makes this type somewhat similar to the movable sleeve type. Not recommended for samples of high viscosity, as the sleeve cannot be cleaned.
Double junction	Combination of the ceramic type and the movable sleeve type overcomes the disadvantages of using either separately. When the outflow of the KCl in the internal solution presents a problem, placing the sample or other salt solution in the external tube will ensure stable measurements.

pH ELECTRODES (3-in-1 ELECTRODES)

Combination electrodes are a glass electrode and a reference electrode incorporated into one unit. 3-in-1 electrodes incorporate a glass electrode and a reference electrode-plus a temperature compensation electrode-into a single unit.

These electrodes are compact and easy to use; they give superb results in pH measurements over a broad range of sample liquids and test conditions. Also, since the glass membrane and the liquid junction are adjacent, only a small amount of sample fluid is required and they are extremely simple to clean. The internal reference electrode uses a solution of 3.33 mol/L KCl.







3-in-1 ELECTRODES **ToupH**

<ToupH> electrode: HORIBA's glass membrane molding technology achieves strengths more than 10 times the Japanese Industrial Standards (strength tests). New dome-shaped construction (9615-10D) boosts strength in all directions.




Type	Applicable temperature range(°C)	pH range	Liquid junction	Internal solution	Feature
9615-10D General laboratory application ToupH Standard ToupH electrode	0-100	0-14	Ceramic	#300 (KCl)	Quick stability, and reduction of drift. No more worries about the timing of your measurement value readings. ●Uses responsive glass that is 10 times stronger than JIS standards. The domed shape provides strength in all directions, greatly reducing damage concerns. ●Constructed with smooth surfaces for easy wiping and cleaning. ●Waterproof ●Pb free (Recommended) Perfect for preparing buffers. Can be used on a wide range of aqueous test solutions. (Post 9611-10D, 6366-10D model)
9618-10D Precious trace amount sample ToupH Micro ToupH electrode	0-60	0-14	Ceramic	#300 (KCl)	This pH electrode with temperature compensation sensor can take measurements from samples as small as 50 µL. ●Compatible with extremely small containers such as micro tubes etc. ●Waterproof ●The temperature sensor is placed next to the response section for high-speed temperature response. (Recommended) Can be used for a wide range of aqueous solutions, including those that cannot be obtained in large quantities. We recommend using our specialized cleaning solution after measuring samples that contain proteins. (Post 9669-10D model)
9680-10D For large containers and long test tubes ToupH Long ToupH electrode	0-100	0-14	Ceramic	#300 (KCl)	283 mm length & 8 mm diameter. The long, thin design makes this electrode perfect for measuring in large containers and test tubes. ●Uses responsive glass that is 10 times stronger than JIS standards. ●Constructed with smooth surfaces for easy wiping and cleaning. ●Waterproof ●Pb free (Recommended) For measuring samples such as microbe culture fluids in test tubes. We recommend that it be used with the long type electrode stand (FA-70L). (Post 9678-10D, 6378-10D model)
9681-10D High viscosity application ToupH Sleeve ToupH electrode	0-60	0-14	Movable sleeve	#300 (KCl)	Stable measurement can also be achieved for highly viscous samples. ●The liquid junction section is constructed with a moveable sleeve that can be rinsed clean, preventing highly viscous samples from clogging the liquid junction, and maintaining stable measurement performance. ●Waterproof ●Pb free (Recommended) For highly viscous samples and solutions, and samples that contain non-aqueous solvents (such as cosmetics or paints). We recommend that you take measurements while using the graph display function to confirm stable responses. (Post 9677-10D model)

3-in-1 ELECTRODES **Plastic Body**

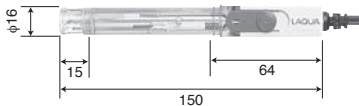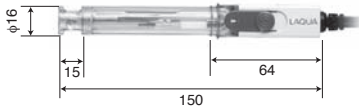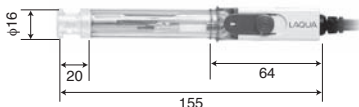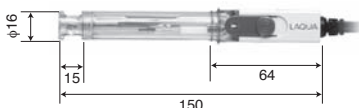Adopting plastic for the body material and covering the electrode tip with a protective tube, this electrode series is ideally suited to measurements in the field and harsh environments. This plastic body lineup consists of the electrodes equipped with high purity glass for tap water measurements, and those with special resistive glasses for hydrofluoric acid and strong alkali sample measurements.



Type	Applicable temperature range(°C)	pH range	Liquid junction	Internal solution	Feature
9625-10D Standard type  3200360505	0-100	0-14	Ceramic	#300 (KCl)	Cased in a plastic body to enable field measurements. The slide-type internal solution filler permits submerged measurements in depths up to 1m (for up to 30 minutes) ● Waterproof ● Pb free (Recommended) Suitable for measurements for tap water, drinking water, field measurements. (Post 9621-10D model)
9630-10D For tap water  3200528726	0-100	0-14	Ceramic	#300 (KCl)	Can measure samples with low conductivity or buffering capacity such as tap water, by adopting high purity multicomponent lithium series glass to its body. Optimal for quality control in Water purification plant. ● Waterproof ● Pb free * Recommended to use along with the dedicated wash solution (model: 230)
9631-10D Hydrofluoric acid resistant model  3200524119	0-60	2-12	Ceramic	#300 (KCl)	Long life capable of measuring about 1000 times*. Rolled glass architecture achieves easy maintenance and long-term reliable measurement. Compliant with the Measurement Act Certification(Japan). Optimal for drain water control after etching process etc. ● Waterproof ● Pb free * When a measurement is conducted for 1 minute with 1% hydrofluoric acid solution (at 25°C)
9632-10D Strong alkali resistant model  3200524120	0-100	0-14	Ceramic	#300 (KCl)	An alkali resistant glass membrane achieves higher resistance and about five times* longer stability than our conventional products. Suitable for strong alkali samples such as plating solutions. ● Waterproof ● Pb free * With 0.1 mol/L sodium (about pH 13) (at 60°C)

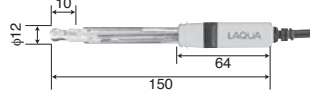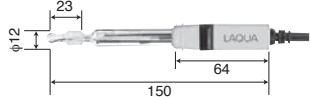3-in-1 ELECTRODES **Glass Body**

Type	Applicable temperature range(°C)	pH range	Liquid junction	Internal solution	Feature
6367-10D Standard type (sleeve)  3014079136 (9003011800)	0-60	0-14	Sleeve	#300 (KCl)	Uses a sleeve for the liquid junction, improving the stability and repeatability. For measuring pH at high accuracy. (Standard accessory for model F-24II.)
6377-10D For measurement of low-conductivity water and non-aqueous solvents  3014093085 (9003014100)	0-60	0-14	Movable sleeve	#300 (KCl)	Uses a glass membrane highly sensitive to low-conductivity water and non-aqueous solvents. Movable sleeve used at the liquid junction.
6252-10D For food application (needle type)  3014080850 (9003013800)	0-60	0-12	Ceramic	#300 (KCl)	Needle electrode allows measurement of aqueous solutions too.

pH ELECTRODES (ISFET (Semiconductor electrode))

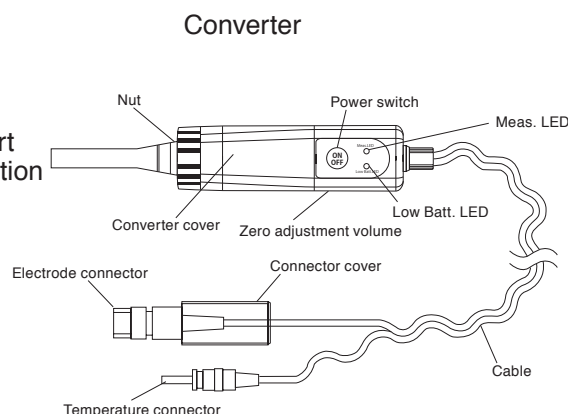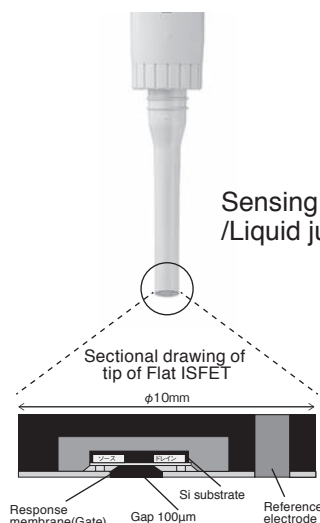ISFET is the abbreviation of Ion Sensitive Field Effect Transistor. Since ISFET is robust and will not crack like the conventional glass electrodes, it can be easily handled and maintained. The response part is equipped with a flat and miniature semiconductor-based sensor, which makes the measurement even on extremely small samples possible. Combination of HORIBA's unique semiconductor device structure and improvement of the electrostatic protection circuit enables to reduce greatly the static electricity effect that had been the weak point of the semiconductor sensor. Now the measurement has become more comfortable and reliable.



ISFET ELECTRODES ISFET

Type	Applicable temperature range(°C)	pH range	Liquid junction	Feature
0040-10D Surface of solid samples Flat ISFET pH electrode 0040-10D 	0-60	0-14	Porous sintered polyethylene	<p>The sensor is located on the flat surface of the tip, with less than a 100 µm difference from the housing.</p> <ul style="list-style-type: none"> ●Measurements can be made from a minute amount of moisture on the solid sample surface. ●Use of a semiconductor sensor means there are no concerns that the electrode will be damaged. ●Also perfect for measuring samples in shallow containers such as Petri dishes. ●Waterproof ●Replaceable Sensor (Recommended) <p>For surface measurement of gelatinous materials such as nutrient agar, and foodstuffs such as meat.</p> <p>Evaluation of sheet materials such as cloth or paper. If the sample only has a small amount of moisture, pure water etc. is required.</p>
0141 Replacement sensor for 0040-10D 3200367925				
0030-10D Inside solid samples Needle ISFET electrode 	0-60	0-14	ABS, epoxy, polyethylene, Ta ₂ O ₅ , platinum	<p>The sharp tip can pierce solid samples to take measurements.</p> <ul style="list-style-type: none"> ●Use of a semiconductor sensor means there are no concerns that the electrode will be damaged. ●Waterproof (Recommended) <p>For measuring inside foodstuffs, such as fruits, vegetables and bread.</p>
0131 Replacement sensor for 0030-10D 3014028323 (9096002100)				

REPLACEMENT SENSOR for ISFET ELECTRODES

Type	Feature
141 Flat ISFET (0040-10D) sensor part (for replacement) 	Replaceable sensor tip for flat ISFET electrode (0040-10D)
131 Needle ISFET (0030-10D) sensor part (for replacement) 	Replaceable sensor tip for needle ISFET electrode (0030-10D)

pH ELECTRODES (GLASS ELECTRODES(G), REFERENCE ELECTRODES(R))

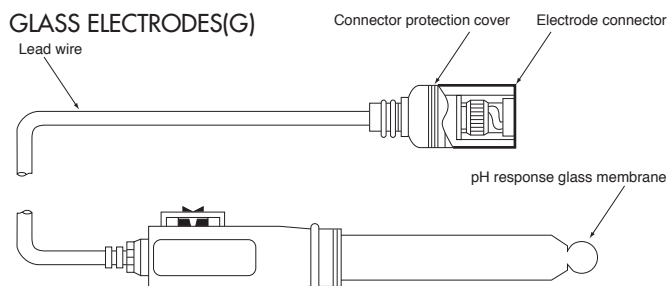
Glass electrodes measure the pH value in the sample solution by detection of electromotive force, i.e., voltage.

HORIBA's superior glass electrodes have all the qualities required for accurate measurement and testing: they are responsive to changes in electromotive force, sensitive to very slight alkaline differences, have a low internal resistance, and are extremely durable. HORIBA's electrodes are perfect not only for laboratory pH measurement conditions, but are in widespread general use for pH measurement.


Our series of electrodes for use with HORIBA's F, M, & D Series of pH meters incorporate a composite lithium glass for the pH-responsive glass membrane. This gives them extremely high sensitivity. They connect to the industry-standard universal BNC connectors. The holder portion has a squared-off design to prevent the electrode from rolling, protecting it from damage.

Reference electrodes constitute part of the detection portion of pH meters; they are used together with a glass electrode to isolate the electromotive force generated in the glass electrode. HORIBA's reference electrodes use a top-quality internal reference electrode and a liquid junction with numerous special features; this gives them an incredible stable indication of electrical potential, making them particularly suitable as reference electrodes in all types of pH and electrical potential measurement. These electrodes have a double-junction configuration, incorporating two types of liquid junction, using capillary tubes, a sleeve with large surface area, and an easy-to-use ceramic filter.

GLASS ELECTRODES(G)



REFERENCE ELECTRODES(R)

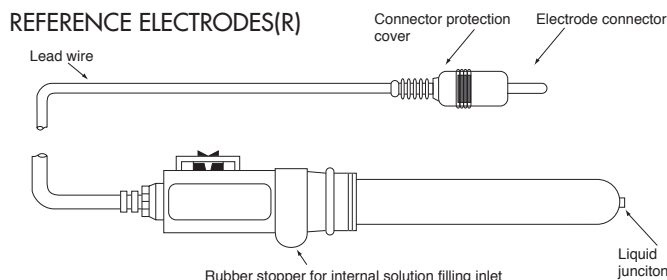

Glass Electrodes(G)



Type	Usage	Applicable temperature range(°C)	pH range	Applicable reference electrode	Feature
1066A-10C Standard type 3014080432 (9003012200)	 Glass electrode 1066A-10C Reference electrode 2060A-10T or other	0-100	0-14	2060A 2565A	Very durable minimum alkali errors. Most widely used for general pH measurements.
1076A-10C For measurement of low-conductivity water and non-aqueous solvents. 3014093084 (9003014200)	 Glass electrode 1076A-10C Reference electrode 2060A-10T or other	0-100	0-14	2060A 2565A	Uses a glass membrane highly sensitive to low-conductivity water and non-aqueous solvents. Can also be used for ordinary pH measurement.

Reference Electrodes(R)


Type	Applicable temperature range(°C)	Liquid junction	Internal solution	Applicable glass electrode	Feature
2060A-10T Standard type 3014080434 (9003012500)	0-100	Ceramic	#300 (KCl)	1066A 1076A	Suitable for a wide range of pH measurements since the resistance of the liquid junction is small.
2565A-10T Double-junction type 3014080436 (9003012700)	0-100	Intermediate: Ceramic External: Sleeve	#300 (KCl)	1066A 1076A	Suitable for measurements of liquid other than normal aqueous solutions, such as suspensions, emulsions, paste, and non-aqueous solutions. When the potassium chloride solution of the internal solution reacts with the sample, measurements can be stably carried out by filling the sample or any other chloride solution in the external jacket. The replacement of the internal solution and the cleaning of the liquid junction can be carried out easily.

pH ELECTRODES (COMBINATION), Temperature Compensation Electrode, METALLIC ELECTRODES (FOR ORP MEASUREMENT)

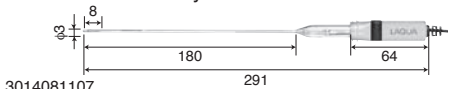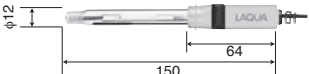
Combination Electrodes

Type	Applicable temperature range(°C)	pH range	Liquid junction	Internal solution	Feature
6069-10C For very slender test tubes  3014081107 (9003013500)	0-60	0-14	Ceramic	#310 (KCl with AgCl)	For measuring pH of a small amount of sample in a slender tube (more than 3.5 mm dia.) such as a NMR test tube.
6261-10C Flat type  3014081807 (9003013700)	0-50	0-12	Sleeve	#300 (KCl)	Since the pH response membrane and the liquid junction are located on the same surface, pH values on the surfaces of skin, leather, paper, and leaves can be measured.

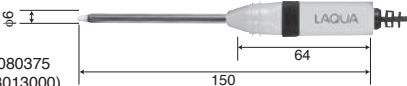Temperature Compensation Electrode

Type	Applicable temperature range(°C)	Applicable	Temperature compensation element	Feature
4163-10T  3014080375 (9003013000)	0-100	Temperature compensation and measurement	Thermistor	Used to automatically compensate the changes in the electromotive force of the pH electrode due to temperatures and also to measure temperatures.

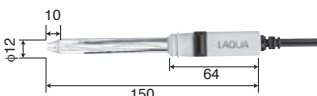Metallic Electrode (For ORP Measurement)


Type	Applicable temperature range(°C)	Electrode material	Applicable reference electrode	Internal solution	Feature
9300-10D Waterproof platinum combination type  3014046710 (9096000400)	0-60	Pt	—	#300 (KCl)	Waterproof. Uses a flat type metallic electrode, which allows a small amount of sample to be measured.

ION ELECTRODES

See P4 for information about reference electrodes.











Ion-selective electrodes are responsive to concentration of particular ions in the test liquid and are variable-potential electrodes. They are used in conjunction with reference electrodes to measure the concentration of particular ions. HORIBAs years of experience and know-how in this field are behind the wide range of ion electrodes we offer.



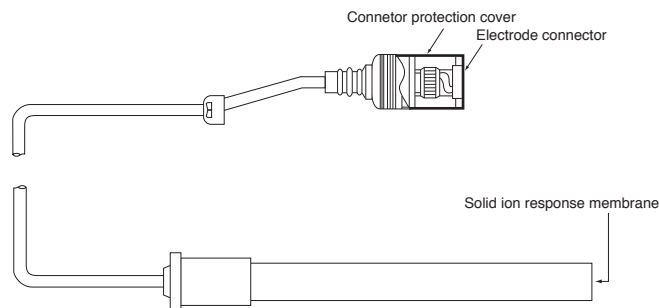
When measurements are made using an ion meter, by calibrating with various standard solutions, direct readings of the concentration of the ion in question can be taken. Note that since volume-detection level changes with temperature, measurements must be taken at a fixed temperature.



①: Measuring range ②: pH range ③: Applicable temperature range ④: Response time (90%)

Type	Measuring range	Applicable reference electrode	Selection coefficient
Cyanide ion electrode 8001-10C 3014094393 (9003015500)	①: 0.03 to 2,600 mg/L CN ⁻ (10 ⁻⁶ to 10 ⁻¹ mol/L CN ⁻) ②: 2.6 mg/L (10 ⁻⁴ mol/L) CN ⁻ pH 12 to 13 ③: 0 to 50°C ④: Within 10 seconds	2060A, 2565A	S ²⁻ , MnO ₄ ⁻ = Not acceptable I ⁻ = 0.1 S ₂ O ₃ ²⁻ = 1
Chloride ion electrode (combination) 6560-10C 3014093430 (9003014500)	①: 0.35 to 35,000 mg/L Cl ⁻ (10 ⁻⁵ to 1 mol/L Cl ⁻) ②: 350 mg/L (10 ⁻² mol/L) Cl ⁻ pH 3 to 11 ③: 0 to 50°C ④: Within 5 seconds	—	S ₂ O ₃ ²⁻ , S ²⁻ , I ⁻ , Ag ⁺ , Hg ²⁺ = Not acceptable SCN ⁻ = 0.3 MnO ₄ ⁻ = 0.1 Br ⁻ = 0.03 NO ₃ ⁻ , F ⁻ , HCO ₃ ⁻ , SO ₄ ²⁻ , PO ₄ ²⁻ = 1,000
Chloride ion electrode 8002-10C 3014094394 (9003015600)	①: 0.35 to 35,000 mg/L Cl ⁻ (10 ⁻⁵ to 1 mol/L Cl ⁻) ②: 350 mg/L (10 ⁻² mol/L) Cl ⁻ pH 3 to 11 ③: 0 to 50°C ④: Within 5 seconds	2565A	S ₂ O ₃ ²⁻ , S ²⁻ , I ⁻ , Ag ⁺ , Hg ²⁺ = Not acceptable SCN ⁻ = 0.3 MnO ₄ ⁻ = 0.1 Br ⁻ = 0.03 NO ₃ ⁻ , F ⁻ , HCO ₃ ⁻ , SO ₄ ²⁻ , PO ₄ ²⁻ = 1,000
Sulfide ion electrode 8003-10C 3014094395 (9003015700)	①: 0.32 to 32,000 mg/L S ²⁻ (10 ⁻⁵ to 1 mol/L S ²⁻) ②: 3.2 mg/L (10 ⁻⁴ mol/L) S ²⁻ pH 12 to 14 ③: 0 to 50°C ④: Within 10 seconds	2060A, 2565A	CN ⁻ = Not acceptable S ₂ O ₃ ²⁻ = 10 I ⁻ , F ⁻ , Cl ⁻ , PO ₄ ²⁻ , SO ₄ ²⁻ = 1,000
Iodide ion electrode 8004-10C 3014094396 (9003015800)	①: 0.0127 to 12,700 mg/L I ⁻ (10 ⁻⁷ to 10 ⁻¹ mol/L I ⁻) ②: 1,270 mg/L (10 ⁻² mol/L) I ⁻ pH 2 to 11 ③: 0 to 50°C ④: Within 10 seconds	2060A, 2565A	MnO ₄ ⁻ , S ²⁻ , CN ⁻ = Not acceptable S ₂ O ₃ ²⁻ = 10 NO ₃ ⁻ = 100 Br ⁻ = 1,000
Bromide ion electrode 8005-10C 3014094397 (9003015900)	①: 0.8 to 80,000 mg/L Br ⁻ (10 ⁻⁵ to 1 mol/L Br ⁻) ②: 800 mg/L (10 ⁻² mol/L) Br ⁻ pH 1.5 to 11.5 ③: 0 to 50°C ④: Within 5 seconds	2565A	S ₂ O ₃ ²⁻ , I ⁻ , S ²⁻ , CN ⁻ = Not acceptable MnO ₄ ⁻ = 1 Cl ⁻ , PO ₄ ²⁻ = 100 F ⁻ , NO ₃ ⁻ , SO ₄ ²⁻ = 1,000
Copper ion electrode 8006-10C 3014094398 (9003016000)	①: 0.06 to 6,350 mg/L Cu ²⁺ (10 ⁻⁶ to 10 ⁻¹ mol/L Cu ²⁺) ②: 6.35 mg/L (10 ⁻⁴ mol/L) Cu ²⁺ pH 2 to 6 ③: 0 to 50°C ④: Within 10 seconds	2565A	Fe ²⁺ = 0.1 Ni ²⁺ , Na ⁺ = 1,000
Cadmium ion electrode 8007-10C 3014094399 (9003016100)	①: 0.1 to 11,240 mg/L Cd ²⁺ (10 ⁻⁶ to 10 ⁻¹ mol/L Cd ²⁺) ②: 11 mg/L (10 ⁻⁴ mol/L) Cd ²⁺ pH 3 to 8 ③: 0 to 50°C ④: Within 10 seconds	2060A, 2565A	Cu ²⁺ , Hg ²⁺ , Ag ⁺ = Not acceptable Pb ²⁺ = 0.1 Fe ³⁺ = 1 Cr ³⁺ , Fe ²⁺ = 100 Ni ²⁺ = 1,000

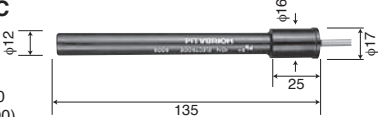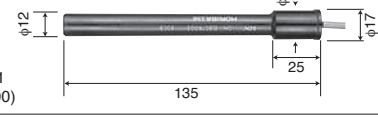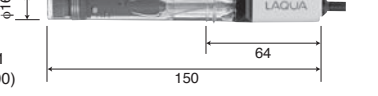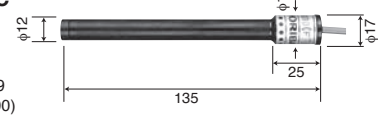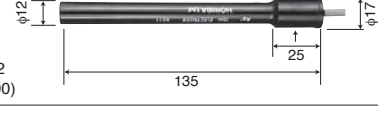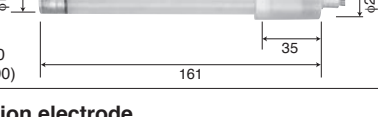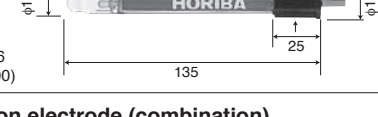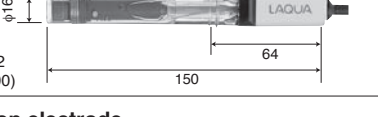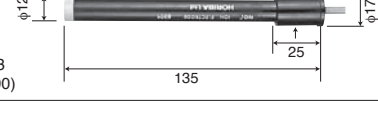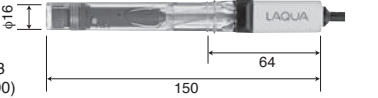Type	Measuring range	Applicable reference electrode	Selection coefficient
Lead ion electrode 8008-10C  3014094400 (9003016200)	①: 2 to 20,000 mg/L Pb ²⁺ (10 ⁻⁵ to 10 ⁻¹ mol/L Pb ²⁺) ②: 20 mg/L (10 ⁻⁴ mol/L) Pb ²⁺ pH 4.5 to 6.5 ③: 0 to 50°C ④: Within 10 seconds	2565A	Cu ²⁺ , Hg ²⁺ , S ²⁻ , Ag ⁺ = Not acceptable Fe ³⁺ = 0.01 Cr ³⁺ = 1 Cd ²⁺ = 10 Ni ²⁺ , Mg ²⁺ , Zn ²⁺ = 100 NH ₄ ⁺ , K ⁺ = 1,000
Thiocyanate ion electrode 8009-10C  3014094401 (9003016300)	①: 0.6 to 5,800 mg/L SCN ⁻ (10 ⁻⁵ to 10 ⁻¹ mol/L SCN ⁻) ②: 5.8 mg/L (10 ⁻⁴ mol/L) SCN ⁻ pH 2 to 12 ③: 0 to 50°C ④: Within 30 seconds	2565A	CN ⁻ , I ⁻ , S ²⁻ , S ₂ O ₃ ²⁻ = Not acceptable Br ⁻ = 1 Cl ⁻ = 100
Fluoride ion electrode (combination) 6561-10C  3014093431 (9003014600)	①: 0.02 to 19,000 mg/L F ⁻ (10 ⁻⁶ to 1 mol/L F ⁻) ②: 20 mg/L (10 ⁻³ mol/L) F ⁻ pH 4 to 10 ③: 0 to 50°C ④: Within 5 seconds	—	Possible interference when multiply-charged ion (ex. Al ³⁺ , Fe ³⁺) coexisted and foamed the complex.
Fluoride ion electrode 8010-10C  3014093439 (9003016400)	①: 0.02 to 19,000 mg/L F ⁻ (10 ⁻⁶ to 1 mol/L F ⁻) ②: 20 mg/L (10 ⁻³ mol/L) F ⁻ pH 4 to 10 ③: 0 to 50°C ④: Within 5 seconds k1	2060A, 2565A	Possible interference when multiply-charged ion (ex. Al ³⁺ , Fe ³⁺) coexisted and foamed the complex.
Silver ion electrode 8011-10C  3014094402 (9003016500)	①: 0.01 to 110,000 mg/L Ag ⁺ (10 ⁻⁷ to 1 mol/L Ag ⁺) ②: 1 mg/L (10 ⁻⁵ mol/L) Ag ⁺ pH 2 to 10 ③: 0 to 50°C ④: Within 10 seconds	2565A	Hg ²⁺ = Not acceptable Cu ²⁺ , Cd ²⁺ , Pb ²⁺ , Zn ²⁺ , Mg ²⁺ , Ca ²⁺ , Na ²⁺ , K ⁺ = over 1,000
Ammonia ion electrode (combination) 5002A-10C  3014093560 (9003016600)	①: 0.1 to 1,000 mg/L NH ₃ ②: Adjust more than pH 12 ③: 0 to 50°C ④: Within 30 seconds when substituting low concentration to high concentration Within 2 minutes when substituting high concentration to low concentration	—	—
Sodium ion electrode 1512A-10C  3014068526 (9003016700)	①: 2.3 to 230,000 mg/L Na ⁺ (10 ⁻⁴ to 10 mol/L Na ⁺) ②: 230 mg/L (10 ⁻² mol/L) Na ⁺ Over pH 4.5 ③: 0 to 60°C ④: Within 30 seconds k1	2565A	K ⁺ , Li ⁺ = 10 NH ₄ ⁺ = 20 Ca ²⁺ = 500
Nitrate ion electrode (combination) 6581-10C  3014093432 (9003014700)	①: 0.62 to 62,000 mg/L NO ₃ ⁻ (10 ⁻⁵ to 1 mol/L NO ₃ ⁻) ②: 62 mg/L (10 ⁻³ mol/L) NO ₃ ⁻ pH 3 to 7 ③: 0 to 50°C ④: Within 15 seconds k2	—	ClO ₄ ⁻ = 0.03 I ⁻ = 0.1 Br ⁻ = 2 NO ₂ ⁻ = 3 Cl ⁻ = 40 F ⁻ = 200 CH ₃ COO ⁻ = 300 SO ₄ ²⁻ = over 1,000
Nitrate ion electrode 8201-10C  3014094403 (9003016800)	①: 0.62 to 62,000 mg/L NO ₃ ⁻ (10 ⁻⁵ to 1 mol/L NO ₃ ⁻) ②: 62 mg/L (10 ⁻³ mol/L) NO ₃ ⁻ pH 3 to 7 ③: 0 to 50°C ④: Within 15 seconds k2	2565A	ClO ₄ ⁻ = 0.03 I ⁻ = 0.1 Br ⁻ = 2 NO ₂ ⁻ = 3 Cl ⁻ = 40 F ⁻ = 200 CH ₃ COO ⁻ = 300 SO ₄ ²⁻ = over 1,000
Potassium ion electrode (combination) 6582-10C  3014093433 (9003014800)	①: 0.04 to 39,000 mg/L K ⁺ (10 ⁻⁶ to 1 mol/L K ⁺) ②: 3.9 mg/L (10 ⁻⁴ mol/L) K ⁺ pH 5 to 11 ③: 0 to 50°C ④: Within 15 seconds *3	—	Rb ⁺ = 0.4 Cs ⁺ = 3 NH ₄ ⁺ = 70 Li ⁺ , Na ⁺ , Mg ²⁺ , Ca ²⁺ , Sr ²⁺ , Ba ²⁺ = over 1,000

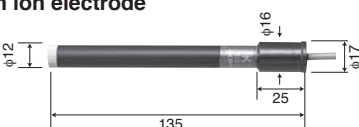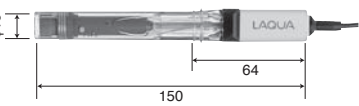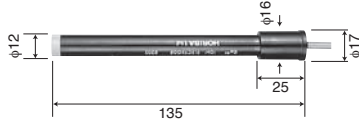Type	Measuring range	Applicable reference electrode	Selection coefficient
Potassium ion electrode 8202-10C  3014094404 (9003016900)	①: 0.04 to 39,000 mg/L K ⁺ (10 ⁻⁶ to 1 mol/L K ⁺) ②: 3.9 mg/L (10 ⁻⁴ mol/L) K ⁺ pH 5 to 11 ③: 0 to 50°C ④: Within 15 seconds k3	2565A	Rb ⁺ = 0.4 Cs ⁺ = 3 NH ₄ ⁺ = 70 Li ⁺ , Na ⁺ , Mg ²⁺ , Ca ²⁺ , Sr ²⁺ , Ba ²⁺ = over 1,000
Calcium ion electrode (combination) 6583-10C  3014093434 (9003014900)	①: 0.4 to 40,080 mg/L Ca ²⁺ (10 ⁻⁵ to 1 mol/L Ca ²⁺) ②: 4.0 mg/L (10 ⁻⁴ mol/L) Ca ²⁺ pH 5 to 11 ③: 0 to 50°C ④: Within 15 seconds k4	—	Fe ³⁺ = 0.1 Fe ²⁺ , Zn ²⁺ = 1 Sr ²⁺ = 50 Ni ²⁺ , Cu ²⁺ = 70 Co ²⁺ = 350 Mn ²⁺ = 500 Mg ²⁺ = 1,000 Na ⁺ , K ⁺ , Ba ²⁺ , NH ₄ ⁺ = over 1,000
Calcium ion electrode 8203-10C  3014068839 (9003017000)	①: 0.4 to 40,080 mg/L Ca ²⁺ (10 ⁻⁵ to 1 mol/L Ca ²⁺) ②: 4.0 mg/L (10 ⁻⁴ mol/L) Ca ²⁺ pH 5 to 11 ③: 0 to 50°C ④: Within 15 seconds k4	2060A, 2565A	Fe ³⁺ = 0.1 Fe ²⁺ , Zn ²⁺ = 1 Sr ²⁺ = 50 Ni ²⁺ , Cu ²⁺ = 70 Co ²⁺ = 350 Mn ²⁺ = 500 Mg ²⁺ = 1,000 Na ⁺ , K ⁺ , Ba ²⁺ , NH ₄ ⁺ = over 1,000

• Sensor holder is necessary for ion electrode except of combination type to attach to electrode stand.

• The response time is the time which is required to reach 90% response when the ion concentration is gradually changed from 10⁻⁴ mol/L to 10⁻² mol/L with the solution stirred.

Exception:

- *1: 90% response when ion concentration is changed to 10⁻⁶ mol/L ~ 10⁻² mol/L
- *2: 95% response when ion concentration is changed to 10⁻³ mol/L ~ 10⁻¹ mol/L
- *3: 95% response when ion concentration is changed to 10⁻⁴ mol/L ~ 10⁻² mol/L
- *4: 95% response when ion concentration is changed to 10⁻⁴ mol/L ~ 10⁻¹ mol/L


• The selection coefficient is a ratio of the limit concentration of coexisting ions (mol/L) to the ion concentration to be measured (mol/L); The value of 1000 means that the coexisting ions can be permitted up to 1000 times the ion measured and "not acceptable" means that chemical change occurs in the solid response membrane.

Cartridges for Ion Sensor

Type	Feature
7660 Chloride ion cartridge 3014093436(9003015000)	Replacement electrode tip for combination or single electrodes
7661 Fluoride ion cartridge 3014093438(9003015100)	
7681 Nitrate ion cartridge 3014068364(9003015200)	
7682 Potassium ion cartridge 3014069795(9003015300)	
7683 Calcium ion cartridge 3014068795(9003015400)	
Membrane(NH₃) 3014067083(9012001000)	Membrane set (6 pcs) for NH ₃ electrodes
370 Internal solution for NH₃ electrodes 3014067184(9012000900)	Contains 250 mL
O-ring 3200043723(9012001100)	Neoprene ring set (10 pcs) for NH ₃ electrodes (JIS B 2401-P7)

CONDUCTIVITY ELECTRODE CELLS

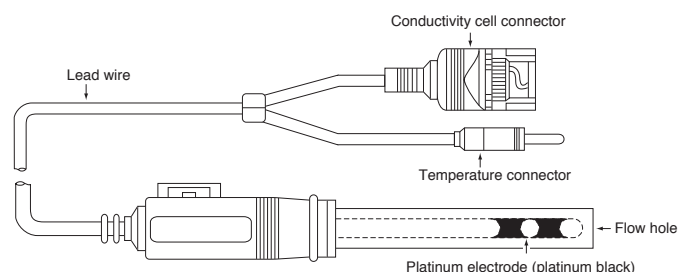Conductivity is calculated as the inverse of the resistance R(in ohms) of the sample solution as $S/m = V/m$ between two parallel electrode plates with a surface area of $1m^2$ separated by a distance of $1m$. Since conductivity changes depending on temperature of the sample solution, values are shown at the standard temperature equivalent of $25^{\circ}C$. HORIBA's conductivity electrodes also have a built-in thermistor for temperature measurement, making them perfect for temperature measurement and for obtaining values equivalent to those at the standard $25C$, when used in conjunction with the conductivity meter. Since the conductivity gives valuable information about the ion composition of the sample solution, it is expect that these useful electrodes will continue to find a wide range of applications in the future.



Conductivity Cells (Submersible Type)

(*1) The cell constants are within 10% of the values shown.

Type	Cell constant (cm ⁻¹)	Measuring range	Sample amount required (mL)	Temperature compensation element	Applicable temperature range (°C)	Remarks
3551-10D 3014081712 (9056000800)	0.1	10 μS/m to 1 S/m (0.1 μS/cm ~10 mS/cm)	50	Incorporated	0-60	For low conductivity water (deionized water or other)
3552-10D 3014081545 (9056000900)	1	0.1 mS/m to 10 S/m (1 μS/cm ~100 mS/cm)	15	Incorporated	0-100	For general purposes (provided as a standard accessory for the DS-10 series)
3553-10D 3014081714 (9056001000)	10	1 mS/m to 100 S/m (10 μS/cm ~1 S/cm)	50	Incorporated	0-60	For high conductivity water
9382-10D 3014046709 (9096000300)	1	0.1 mS/m to 10 S/m (1 μS/cm ~100 mS/cm)	20-30	Incorporated	0-80	Waterproof. For general purposes.


Conductivity Cells (Flow Type)

(*1) The cell constants are within 10% of the values shown.




Type	Cell constant (cm ⁻¹)	Measuring range	Sample amount required (mL)	Temperature compensation element	Applicable temperature range (°C)	Remarks
3561-10D 3014082350 (9056001100)	0.1	10 mS/m to 1 S/m (0.1 μS/cm ~10 mS/cm)	10	Incorporated	0-60	For low conductivity water (pure water or other)
3562-10D 3014082513 (9056001200)	1	0.1 mS/m to 10 S/m (1 μS/cm ~100 mS/cm)	16	Incorporated	0-60	For general purposes
3573-10C 3014082590 (9056001300)	10	1 mS/m to 100 S/m (10 μS/cm ~1 S/cm)	4	Not provided	0-60	For high conductivity water
3574-10C 3014082592 (9056001400)	10	1 mS/m to 10 S/m (10 μS/cm ~100 mS/cm)	0.25	Not provided	0-80	For column chromatography using a very small amount of sample

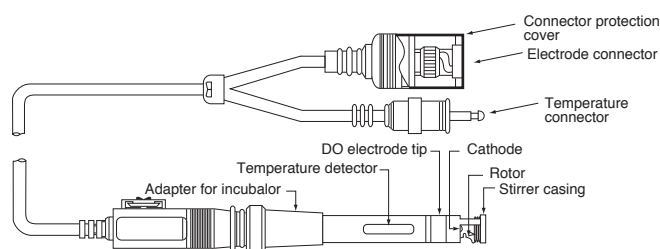DISSOLVED OXYGEN(DO) ELECTRODE & TIPS

Dissolved Oxygen(DO) electrode detect oxygen that diffuses through the oxygen-permeable membrane to determine the amount of dissolved oxygen. The method for measuring dissolved oxygen based on this principle is referred to as the diaphragm electrode method. DO measurement can be carried out much more simply than chemical analysis, which requires complex preparatory procedures to eliminate the effects of deoxidized and oxidized substances. HORIBA's DO electrodes use innovative disposable probe tips. This eliminates the troublesome replacement of membranes and fluid that plagued conventional methods. Each disposable tip comes with its own rotor, so it is not necessary to prepare a separate rotor for each sample. In addition, the electrode has an adaptor for easy use with an incubator in BOD measurement.

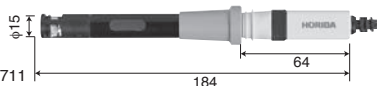
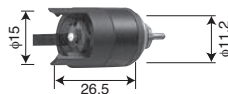


Dissolved Oxygen Electrodes

Type	Applicable temperature range(°C)	Measuring range	Response time	Feature
9520-10D For laboratories  3014046711 (9096000500)	0-45	DO: 0-19.99mg/L Temperature: 0-40°C (When used with dissolved D-25)	20 seconds (90% response time at constant temperature)	Waterproof. Uses a thermistor with a disposable ship-type electrode 7541 as the thermometric element.
9551-20D For field immersible type (2 m cable)  3014047090 (9096002300)	0-40	DO: 0-19.99mg/L Temperature: 0-40°C (When used with dissolved D-55, OM-51)	30 seconds (90% response time at constant temperature)	Waterproof. Uses a thermistor with a disposable ship-type electrode 5401 as the thermometric element.
9551-100D For field immersible type (10 m cable)  3014047091 (9096002400)	0-40	DO: 0-19.99mg/L Temperature: 0-40°C (When used with dissolved D-55, OM-51)	30 seconds (90% response time at constant temperature)	Waterproof. Uses a thermistor with a disposable ship-type electrode 5401 as the thermometric element.







Dissolved Oxygen Electrode Tips

*A commercially available stirrer should be used.

Type	Remarks
5401  3014072770 (9033010000)	A DO electrode chip for replacement. (For the above-mentioned 9551-20D, 9551-100D, 9550-20D, 9550-100D, 5450-20D and 5450-100D)
7541  3014074145 (9074000200)	A DO electrode chip for replacement. (For the above-mentioned 5410-10C, 9520-10D)

ACCESSORIES

For Electrode

Sensor holder (2 pcs/pack) 3200373961  For attaching an ion electrode or the like with a round electrode cap to the stand arm.	9621, 9625 Electrode protector tube (5 pcs/pack) 3200044409 (9003012000)  Protects the tip of the 9621-10D, 9625-10D electrode. Because the electrode is already encased in a plastic sheath, just slip this protective tube over the tip and your pH meter is ready to for work in the field.	Electrode protector cap (5 pcs/pack) 3200043508 (9003012100)  Protects electrodes during storage or transportation.
Electrode protector cap (3 pcs/pack) 3200382477  Protects electrodes during storage or transportation for 9615-10D, 9618-10D, 9681-10D.	Electrode protector cap for long electrode 32000382482  Protects electrodes during storage or transportation for 9680-10D and 9678-10D.	Plug for internal solution filler port (3 pcs/pack) 3200382468  For 9615-10D, 9618-10D, 9681-10D, 9680-10D.





Meter Accessories Meter and Accessories table

	Model	Part No.	Meter									
			F-50	F-70	DS-50	DS-70	D-50	D-70	ES-50	ES-70	OM-50	OM-70
Printer	Printer (for GLP/GMP compliance)	3014030145	○	○	○	○	○*1	○*1	○	○	○	○
	Printer cable	3014030148	○	○	○	○	○*1	○*1	○	○	○	○
	Printer paper	3014030149	○	○	○	○	○*1	○*1	○	○	○	○
	Ink ribbon	3014030150	○	○	○	○	○*1	○*1	○	○	○	○
For Inspection	Digital simulator X-51	—	○	○	—	—	○	○	—	—	○	○
	Digital simulator X-52	—	○*2	○*2	○	○	○*2	○*2	○	○	—	—
Communication and Output	USB cable	3200373941	—	○	—	○	—	—	—	—	—	—
	Serial cable	3014030151	○	○	○	○	○*1	○*1	○	○	○	○
	Analog (alarm) output cable	3014030152	○*3	○*3	○*3	○*3	—	—	—	—	—	—
	COMPACTFLASH® memory card	3014030160	○*4	—	—	—	—	—	—	—	—	—
Power	Data Collection Software *5	—	○	○	○	○	○*1	○*1	○	○	○	○
	AC adapter and cable set	—	○	○	○	○	○	○	○	○	○	○
Meter Accessories	LCD protection sheet	3200382462	—	○	—	○	—	—	—	—	—	—
	Protection cover	3200382441	—	○	—	○	—	—	—	—	—	—
	Electrode hook	3200528475	—	—	—	—	—	○	—	○	—	○
Electrode stand	Electrode stand for F-50/DS-50	3014028342	○	—	○	—	—	—	—	—	—	—
	FA-70S Electrode stand (adjustable type)	3200382557	○	○	○	○	○	○	○	○	○	○
	FA-70L Electrode stand (long type)	3200382560	○	○	○	○	○	○	○	○	○	○
	DP-50S Electrode stand	3014028590	—	—	—	—	○	—	○	—	○	—
	DP-70S Electrode stand	3200528474	—	—	—	—	—	○	—	○	—	○



*1 Except D-51/D-71 *2 Conductivity measurement model: F-54/F-55/F-74/F-74BW/D-54/D-74 *3 Except F-51/F-71/F-74BW/DS-51/DS-71

*4 Only F-53/F-54/F-55 *5 Data collection software is available as a free download for registered users. <http://www.horiba.co.jp/register>


Printer-related

Printer (for GLP/GMP compliance)	Printer cable	Printer paper	Ink ribbon
CBM-910-24RJ100-A There are printers for 100V, 120V and 230V power supplies. Please consult our sales staff when ordering 120V and 230V models.  The model numbers for 120V and 230V are listed below. 120V: CBM-910-24RJ-120-A (3014030146) 230V: CBM-910-24RJ-230-A (3014030147)	3014030148 (9096003800) Cable to connect Printer with 50 series and 70 series. 	3014030149 (9096003900) 20 rolls 	3014030150 (9096004000) 5 pcs/set 





For Inspection

Digital simulator X-51	Digital simulator X-52
pH, mV, ION, DO simulator (for periodic inspection of the electrode) 	Conductivity simulator (for periodic inspection of the electrode) 


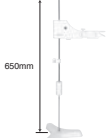





Communication/output

USB cable	Serial cable	COMPACTFLASH® memory card	Analog (alarm) output cable
3200373941  Cable to connect a meter and PC.	3014030151 (9096004800)  Cable to connect a meter and PC (Serial, 9 pins)	3014030160 (9096003000)  For F-53, 54, 55	3014030152 (9096004900)  For F-52, 53, 54, 55 and DS-52, F-72, F-73, F-74, DS-72

Power supply and Meter accessories


AC adaptor cable set.	LCD protection sheet (2 pcs/pack)	Protection cover	Electrode hook
AC adaptor 1.8m cable 1m  120V: 3014031951 230V: 3014031952	3200382462  For F-70, DS-70 series	3200382441  Protects the meter for F-70, DS-70 series	3200528475 For D-70, ES-70, OM-70  With electrode cable winding function *Meter and electrode are not included

Electrode stands


Electrode stand for F-50/DS-50	FA-70L Electrode stand (long type)	FA-70S Electrode stand	DP-50S Electrode stand
3014028342 (9096002600)  For F-50, DS-50	3200382560  Electrode stand (adjustable type)	3200382557  Electrode stand (adjustable type)	3014028590 (9096002700)  For D-50, ES-50, OM-50 *Meter and electrode are not included
DP-70S Electrode stand	Arm, for electrode stand (adjustable type)	Arm, for electrode stand	
3200528474 For D-70, ES-70, OM-70  *Meter and electrode are not included	3200373991  For FA-70S, FA-70L, DP-70S and FA-20S. Also available for FA-50S and "Electrode stand for F-50/DS-50"	3014030158 (9096002800)  For "DP-50S", "FA-50S" and "Electrode stand for F-50/DS-50"	

■ Dimension Unit: mm

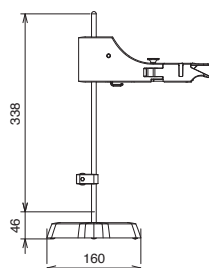
FA-70S




FA-70L

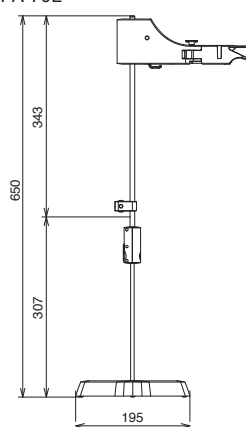

DP-70S

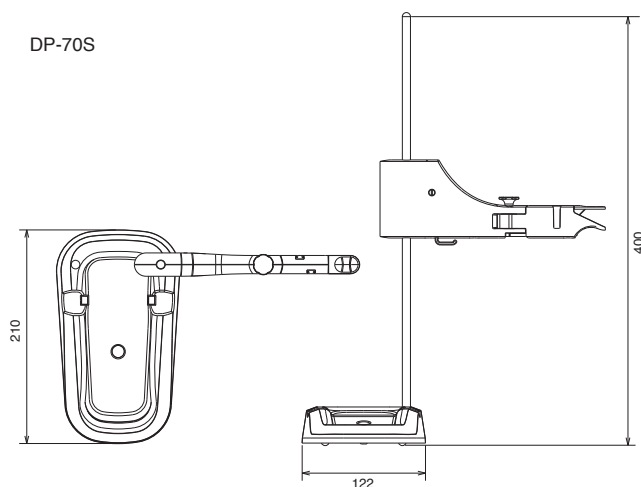


Maintenance Parts for Obsolete Models


Output cord	AC-10 AC adapter	Printer paper (10 rolls)	Dual electrode holder
3200044408 (9078000200) Connect a recorder to make easy work of data analysis after measurement. Applicable models: D-20, 10, OM-10 and D-10 series	3200044196 (9078000100) Applicable models: D-20, F-20, ES-10, OM-10, D-10 and DS-10 series	3200043956 (9079000400) Applicable models: F-15, 16, DS-15, and F-20 series	3200043613 (9096001100)  Applicable model: D-20 series Adaptor for fitting two electrodes

STANDARD SOLUTIONS, INTERNAL SOLUTION for REFERENCE ELECTRODE & CLEANING SOLUTIONS


pH Standard Solution SET (accuracy: ± 0.02 pH)

	Type	Name	pH value(25°C)	Volume(mL)	Remarks
	101-S 3200043642 (9003003500)	Phosphate standard equimolal solution	6.86	500	Use undiluted. The set contains standard and internal solutions, as shown.
		Phthalate standard solution	4.01	250	
		Borate standard solution	9.18	250	
		Internal Solution for Reference Electrode (300)	—	250	


pH Standard Solution (accuracy: ± 0.02 pH)

	Type	Name	pH value(25°C)	Volume(mL)	Remarks
	100-2 3200043639 (9003001500)	Oxalate standard solution	1.68	500	The original solution should be used as it is. For general use as standard solution sets, 101-S (100-4.7.9 and #310 internal solution) are also available.
	100-4 3200043638 (9003001600)	Phthalate standard solution	4.01	500	
	100-7 3200043637 (9003001700)	Phosphate standard equimolal solution	6.86	500	
	100-9 3200043636 (9003001800)	Borate standard solution	9.18	500	
	100-10 3200043635 (9003001900)	Carbonate standard solution	10.02	500	


Condensed pH Standard Solution (accuracy: ± 0.02 pH)

	Type	Name	pH value(25°C)	Volume(mL)	Remarks
	110-4 3200043626 (9003002300)	Condensed phthalate standard solution	4.01	500	Should be diluted when used. The pH values shown are those obtained when the original solution is diluted with pure water at a volume ratio of 1 to 4. For general use.
	110-7 3200043625 (9003002400)	Condensed phosphate standard equimolal solution	6.86	500	
	110-10 3200043624 (9003002500)	Condensed carbonate standard solution	10.02	500	


Standard Solution for Accurate Measurements (N.B.S., accuracy: ± 0.003 pH)

	Type	Name	pH value		Volume(mL)	Remarks
			25°C	37°C		
	100-B4 3200043630 (9003002000)	Phthalate standard solution	4.008	4.030	500	The original solution should be used as it is. This standard solution is for very accurate measurements based on N.B.S. The pH values shown do not necessarily match with those shown in JIS.
	100-B7 3200043631 (9003002100)	Phosphate standard solution	7.413	7.383	500	
	100-B9 3200043632 (9003002200)	Borate standard solution	9.180	9.082	500	

Powder for pH Standard Solution (accuracy: ± 0.05 pH)


	Type	Name	pH value(25°C)	Remarks
	150-4 3200043619 (9003002700)	Powder for phthalate standard solution	4.01	The pH value shown are those obtained when one packet is dissolved in 500 ml of pure water. One packet contains powder for 500 mL. For use in field at factories (10 packets per set)
	150-7 3200043620 (9003002800)	Powder for neutral phosphate standard solution	6.86	
	150-9 3200043621 (9003002900)	Powder for borate standard solution	9.18	

Powder for ORP Standard Solution (accuracy: ± 15 mV)


	Type	Name	ORP value(25°C)	Remarks
	160-51 3200043618 (9003003100)	Powder for ORP standard solution	89 mV (vs. 3.33 mol/L KCl-AgCl)	The ORP values shown are those obtained when one packet is dissolved in 250 mL of pure water. This standard solution should be used immediately after conditioning and can-not be used for 2 hours or more. (10 packets per set)
	160-22 3200043617 (9003003000)	Powder for ORP standard solution	258 mV (vs. 3.33 mol/L KCl-AgCl)	

Note: The pH standard solution by a reliable manufacturer should be selected because they are used as reference for pH measurements. It is recommended for safety not to use the standard liquid which was allowed to stand for long hours after opening its bottle or which was once used.


Internal Solution for Reference Electrode

	Type	Name	Concentration	Volume(mL)	Remarks
	300 3200043640 (9003003200)	For 6327, 6328, F, M, and D-10 series electrodes	3.33 mol/L KCl	250	The original solution should be used as it is. Powder for internal solution (350) is also available for a large amount of internal solution. (The powder is used by dissolving it in pure water.)
	310 3200043622 (9003003300)	For H-7 and old type pH meter electrodes	3.33 mol/L KCl (AgCl, saturation in normal temp.)	250	

Powder for Internal Solution for Reference Electrode

	Type	Remarks
	350 3200043623 (9003003400)	500g. Dissolve in 2L of pure water.

Electrode Cleaning Solution










	Type	Name	Volume(mL)	Remarks
	220 3014028653 (9096002500)	Electrode cleaning solution	50 x 2 pcs	For removing inorganic sample residues from glass electrodes, and for cleaning liquid junctions
	230 3200530494	Electrode cleaning solution	Solution A 30mL(1 bottle) Solution B 100mL(1 bottle)	For 9630-10D (pH electrode for tap water or low conductivity sample)
	250 3200366771	Electrode cleaning solution	400	For removing protein containing sample residues from glass electrodes, and for cleaning liquid junctions.

ACCESSORIES for U-50, U-20XD, U-10, W-20XD SERIES & INTERNAL SOLUTION for REFERENCE ELECTRODE

U-50 Series Accessories




pH sensor 7112 3014057312 (90370048000)	pH sensor <i>Tou</i>pH 7113 3200170923	ORP sensor 7313 3200170920	DO sensor 7543 3200170924	Reference sensor 7210 3200043582 (90370050000)
Reference tip for 7210 3200043587 (90370051000)	Turbidity sensor 7800 for U-52/U-52G 3200172803	Turbidity sensor 7801 for U-53/U-53G 3200172800	Turbidity sensor 7802 for U-54/U-54G 3200318188	DO membrane cap for U-50 series 3200170194
DO inner fluid 306 50 mL 3200170938				

U-20XD/W-20XD Series Sensors

pH sensor 6230 for U-21/22/23, W-22/23 3014050849 (9037005600)	Dissolved oxygen sensor 5460 for U-21/22/23, W-22/23 3014001152 (9037005800)	Chloride ion sensor * 6522 for W-23 3014050860 (9037006000)
		
Fluoride ion sensor * 6530 for W-23 3014050859 (9037006300)	Ammonia sensor 5012 for W-23 3014050864 (9037006200)	pH/ORP sensor 6280 for U-22, W-22/23 3014050850 (9037005700)
		
Nitrate ion sensor * 6531 for W-23 3014050863 (9037005900)	Calcium ion sensor * 6533 for W-23 3014050861 (9037006100)	Potassium ion sensor * 6532 for W-23 3014050862 (9037006400)
		
#5460 membrane replace kit (50 times) 3014050853 (9037007400)		#5012 membrane replace kit (6 units/set) 3014001155 (9037007000)

* Ion selective electrode cartridge included.

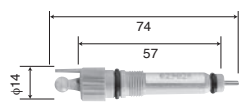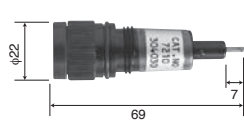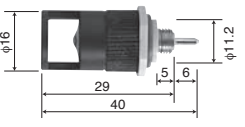U-10 Electrode Tips








pH Electrode Tip 7112 3014057312 (9037004800)		Reference Electrode Tip 7210 3200043582 (9037005000)	
Reference electrode internal solution 330 3200043641 (9037005200) 3.33 mol/L KCL Gel type 250 mL			
Do Tip 7542 3014057313 (9037004900)		Liquid junction 3200043587 (9037005100) For 7210 2 units/set	
Calibration beaker 3200043588 (9037005300)			

Internal Solution for Reference Electrode


300 3200043640 (903003200) Reference electrode internal solution for nitrate ion (250 mL)	301 3014001271 (9037006700) Reference electrode internal solution for chloride (50 mL)	302 3014001273 (9037006600) Reference electrode internal solution for calcium/fluoride (50 mL)	303 3014001272 (9037006900) Reference internal solution for potassium (50 mL)
370 3014067184 (9012000900) Reference internal solution for ammonia (250 mL)	Calibration beaker 3014001156 (9037007300) U-20 For automatic calibration		

ELECTRODES & ACCESSORIES for TWIN/CARDY

SENSOR and ACCESSORIES for LAQUAtwin/TWIN/CARDY

Type	Sample amount required	Measuring temperature	Applicable model	Remarks
S010 LAQUAtwin/TWIN pH sensor  3200459834	Approx. 0.1 mL	5 to 40°C	B-211/B-212 B-213/B-711 B-712/B-713	Liquid junction: Porous macromolecule Glass electrode and reference electrode integrated on a 1mm-thick substrate. Replacement flat type pH sensor.
S021 LAQUAtwin Salt sensor  3200459866	Approx. 0.3 mL	5 to 40°C	B-721	Liquid junction: Porous macromolecule Replacement flat type salt sensor. This sensor respond to sodium ion.
S022 LAQUAtwin Sodium ion sensor  3200459867	Approx. 0.3 mL	5 to 40°C	B-722	Liquid junction: Porous macromolecule Replacement flat type Sodium ion sensor.
S030 LAQUAtwin Potassium ion sensor  3200459868	Approx. 0.3 mL	5 to 40°C	B-731	Liquid junction: Porous macromolecule Replacement flat type Potassium ion sensor.
S040 LAQUAtwin Nitrate ion sensor  3200459870	Approx. 0.3 mL	5 to 40°C	B-341/B-342 B-343/B-741 B-742/B-743	Liquid junction: Porous macromolecule Replacement flat type Nitrate ion sensor.
S050 LAQUAtwin Calcium ion sensor  3200459869	Approx. 0.3 mL	5 to 40°C	B-751	Liquid junction: Porous macromolecule Replacement flat type Calcium ion sensor.
S070 LAQUAtwin Conductivity sensor  3200459672	Approx. 0.12 mL	5 to 40°C	B-771	Replacement flat type Conductivity sensor.






Exclusively for TWIN Conductivity Cell

Type	Measuring range	Cell capacity	Measuring temperature	Temperature compensation element	Remarks
0413 (for B-173)  3014088578 (9088000400)	0 to 19.9mS/cm	Approx. 0.1 mL	5 to 35°C	Incorporated	For B-173 (conductivity meter) only. Cannot be applied for B-771.


















Exclusively for CARDY Ion Electrode

Type	Measuring range	Sample amount required	Measuring temperature	Liquid junction	Remarks
Sodium ion electrode 0221 (for C-121 and C-122) 3014081704 (9076003000)	0.1% (w/w) to 10% (w/w) NaCl	Approx. 0.1 mL	5 to 35°C	Porous macromolecule	For C-121, C-122 (Salt, Sodium ion meter) only. Cannot be applied for B-721, B-722.
Potassium ion electrode 0231 (for C-131) 3014083433 (9076007200)	39 to 3,900 mg/L	Approx. 0.1 mL	5 to 35°C	Porous macromolecule	For C-131 (Potassium ion meter) only. Cannot be applied for B-731.
Nitrate ion electrode 0241 (for C-141) 3014083435 (9076007600)	62 to 6,200 mg/L	Approx. 0.1 mL	5 to 35°C	Porous macromolecule	For C-141 (Nitrate ion meter) only. Cannot be applied for B-340 series & B-740 series.

Accessories

Type	specification	Remarks
Y047 Sampling sheet holder  3200053995	For B-342 (for soil)/replacement sensor (0243) *Cannot be applied for LAQUAtwin B-700 series and their replacement sensors (S010/S021/S022/S030/S040/S050)	For a sample that contain particulate such as soils, suspension. To be used with "Sampling sheet B (model Y046)"
Y048 Sampling sheet holder  3200459736	For LAQUAtwin B-700 series and their replacement sensors (S010/S021/S022/S030/S040/S050) *Cannot be applied for B-342 (for soil)/replacement sensor (0243)	For a sample that contain particulate such as soils, suspension. To be used with "Sampling sheet B (model Y046)"
Y046 Sampling sheet B  3200053858	100 sheets For LAQUAtwin/TWIN series	For trace measurement(0.05 mL),wiping measurement. If a sample that contain particulate, please use *Y047: twin series/Y048: LAQUAtwin series
Y011A Sampling sheet C  3014053435	11 mm x 6 mm x 5 rolls For CARDY series	For trace measurement (0.05 mL), wiping measurement.
Y049 Crop sample press  3200469679	For squeezing a sample such as crop	Standard accessory for B-341, B-741

Standard solution

Type		Value	Volume	Applicable model	Remarks
Y017 Standard solution (pH 6.86) 3200457725		pH 6.86	14 mL 6 bottles	B-711/B-712 B-211/B-212	Replacing model of Y031 (discontinued)
Y014 Standard solution (pH 4.01) 3200457726		pH 4.01	14 mL 6 bottles	B-712/B-212	Replacing model of Y032 (discontinued)
Y021H Standard solution (NaCl 5.0%) 3200457721		NaCl 5.0%	14 mL 6 bottles	B-721/C-121	Replacing model of Y022 (discontinued) *To be used with Y021L for two-point calibration
Y021L Standard solution (NaCl 0.5%) 3200457722		NaCl 0.5%	14 mL 6 bottles	B-721/C-121	Replacing model of Y022 (discontinued) *To be used with Y021H for two-point calibration
Y022H Standard solution (Sodium Ion 2000ppm) 3200457723		Sodium Ion 2000ppm	14 mL 6 bottles	B-722/C-122	Replacing model of Y024 (discontinued) *To be used with Y022L for two-point calibration
Y022L Standard solution (Sodium Ion 150ppm) 3200457724		Sodium Ion 150ppm	14 mL 6 bottles	B-722/C-122	Replacing model of Y024 (discontinued) *To be used with Y022H for two-point calibration
Y031H Standard solution (Potassium Ion 2000ppm) 3200457719		Potassium Ion 2000ppm	14 mL 6 bottles	B-731/C-131	Replacing model of Y025 (discontinued) *To be used with Y031L for two-point calibration
Y031L Standard solution (Potassium Ion 150ppm) 3200457720		Potassium Ion 150ppm	14 mL 6 bottles	B-731/C-131	Replacing model of Y025 (discontinued) *To be used with Y031H for two-point calibration
Y041 Standard solution (Nitrate Ion 5000ppm) 3200053433		Nitrate Ion 5000ppm	14 mL 6 bottles	B-741/B-341	
Y042 Standard solution (Nitrate Ion 300ppm) 3200053514		Nitrate Ion 300ppm	14 mL 6 bottles	B-741/B-742 B-341/B-342	
Y043 Standard solution (Nitrate Ion 2000ppm) 3200053532		Nitrate Ion 2000ppm	14 mL 6 bottles	B-743/B-343 C-141	Replacing model of Y026 (discontinued) *To be used with Y045 for two-point calibration
Y044 Standard solution (Nitrate Ion 30ppm) 3200053535		Nitrate Ion 30ppm	14 mL 6 bottles	B-742/B-342	
Y045 Standard solution (Nitrate Ion 150ppm) 3200053536		Nitrate Ion 150ppm	14 mL 6 bottles	B-743/B-343 C-141	Replacing model of Y026 (discontinued) *To be used with Y043 for two-point calibration
Y051H Standard solution (Calcium Ion 2000ppm) 3200457727		Calcium Ion 2000ppm	14 mL 6 bottles	B-751	
Y051L Standard solution (Calcium Ion 150ppm) 3200457728		Calcium Ion 150ppm	14 mL 6 bottles	B-751	
Y071H Standard solution (Conductivity 12.9mS/cm) 3200457718		Conductivity 12.9mS/cm	14 mL 6 bottles	B-771	
Y071L Standard solution (Conductivity 1.41mS/cm) 3200457717		Conductivity 1.41mS/cm	14 mL 6 bottles	B-771/B-173	Replacing model of Y023 (discontinued) *To be used with Y071H for two-point calibration (Only B-771)

pH Electrode Selection Guide

		3-in-1 ELECTRODES (ToupH)						
		PLASTIC	STANDARD ToupH	LONG ToupH	MICRO ToupH	SLEEVE ToupH	For TAP WATER	HF-PROOF
		9625-10D	9615-10D	9680-10D	9618-10D	9681-10D	9630-10D	9631-10D
Specification	Applicable temperature range (°C)	0-100	0-100	0-100	0-60	0-60	0-100	0-60
	Diameter (mm)	16	12	8	3	12	16	16
	Position of liquid junction (approx.mm)	15	13	21	6	26	15	20
	Length (mm)	150	198	283	185	203	150	155

pH-Sample Conditions

Aqueous Solution	Conductivity	Normal (over 100 mS/m)	●	●	●	●	●	●
		Low (approx.10~100 mS/m)				○	●	
		Very low (approx.5~10 mS/m)				○	○	
		High (approx. 5 S/m)	○	○	○	●	○	○
	Strong alkaline (pH 10-12)			○	○	○		
	Strong acidity (pH 0-2) * Except HF sample			●				●
	Quick heat change (within 50°C)		●				●	●
	High viscosity (approx. 5 Pa·S)					●		
	Containing non-aqueous solvent			○	○	○		
	Suspension			○	○	●		
Solid/Semisolid	Inside							
	Surface							

pH-Sample Conditions

Sample Containers	Microtube/plate (> 50 µL)		×	×	×	●	×	×	×
	NMR tube	φ5 mm ID > φ4 mm	×	×	×	×	×	×	×
	Ampule > φ4 mm					●			
	Micro container (> 2 mL)				○	●			
	Tube	ID:13 mm, L:100 ~ 150 mm			●				
	Beaker	10 mL ~ 1 L	●	●	○	○	○	●	●
	Large container (> 1 L)		○	○	●		○	○	
	Petri dish								
	Droplet		×	×	×	×	×	×	×

pH- Typical Samples

Water	Pure/ion-exchange water (approx. 0.1 mS/m)								
	Distilled water (approx. 0.5 mS/m)			○					
	Tap/drinking water (approx. 10 mS/m)		○	○		○	●		
	Surface water			○		○	●		
	Pharmaceutical water			○		○			
	Enviromental water/acid rain		○	○		○	○		
Chemical reagent/solvent	Caustic/strong acid (Except HF sample)			●		○		●	
	Hydrofluoric acid							●	
	Organic solvent		×				×	×	
	KCl-reactive solution		×	×	×	×	×	×	
	Surfactant			○		●			
	Water-based paint			○		●			
Pharmaceutical/ biology sample	Dye/coloring agent					●			
	Protein-containing sample			○		○	●		
	Medicinal preparation					○			
	Enzyme solution				○	●			
	Tris buffer			●		○	○		
	Suspension			○		●			
Food	Agar medium								
	Jam			○		●			
	Meat/fish								
	Fruit/vegetable								
	Dough								
	Honey								
Beverage/ seasoning	Cheese/butter								
	Yogurt		○	○		○	○		
	Beer		○	○		●	○		
	Milk			○		●			
	Carbonated drink/juice/sauce/soy sauce			○		●			
	Mayonnaise/ketchup			○		●			
Cosmetic/ lotion	Beauty cream/mascara			○		●			
	Gel/soap/shampoo			○		●			
	Hairdye lotion			○		●			
	Emulsified liquid			○		○			

Representative sample names are shown in the table, therefore they may not apply to all cases. A reference electrode is necessary for a glass electrode.

					●				○	○	○	○
					●				○	○	○	○
					○					○		○
					○					○		○
									○		○	○
	×	×	×		○					●		●
	×	×	×	×	×	×	×	×	○	○	×	●
					○				○	○		○
					○				○	○		○
					○				○	○		○
				○						○		○
					○					○		○
						○				○		○
									○	○		○
					●			●	○	○		○
		● (inside)	● (surface)		○	○		○	○	○		○
		● (inside)	● (surface)			●		○				
		● (inside)	● (surface)			○		○				
		● (inside)	● (surface)			○		○				
		○ (inside)	○ (surface)		●					○		○
		● (inside)	○ (surface)			○		○				
		● (inside)	● (surface)	○		○		○	○	○		○
				○	●					●		●
				○	○				○	○		○
				○	○				○	○		○
					○				○	○		○
		○			○	○				○		○
					○					○		○
					○					○		○
					●					○		○

HORIBA WATER QUALITY ANALYZER LINEUP

Benchtop LAQUA

pH/Water Quality Analyzer F-70 series

pH	mV (ORP)	ION	COND	DO	RESI	SAL
TDS						

- Intuitive and very easy to use touch panel operation and navigation
- USB PC Communication and USB memory
- Full support for various country pharmaceutical pure water guidelines (USP/EP/JP/CP)
- Multi-language support (Japanese, English, Chinese, Korean)
- Enhanced data reliability with validation features (GLP/GMP compliance)



DS-70 series

COND	RESI	SAL	TDS
------	------	-----	-----

- Intuitive and very easy to use touch panel operation and navigation
- USB PC Communication and USB memory
- Full support for various country pharmaceutical pure water guidelines (USP/EP/JP/CP)
- Multi-language support (Japanese, English, Chinese, Korean)
- Enhanced data reliability with validation features (GLP/GMP compliance)



Portable LAQUAact

pH/Water Quality Meter D-70 Series

pH	ORP	ION	COND	RESI
SAL	TDS	Oxygen		

- Compact design with special shock and alcohol resistant body
- Fully waterproof/dustproof meter and electrode (IP67 rated)
- Large LCD with backlight
- Data memory up to 1000 sets and PC/Printer output
- Two parameter measurement with simultaneous dual display



Conductivity Meter ES-71

COND	RESI	SAL	TDS
------	------	-----	-----

- Compact design with special shock and alcohol resistant body
- Fully waterproof/dustproof meter and electrode (IP67 rated)
- Large LCD with backlight
- Data memory up to 1000 sets and PC/Printer output
- Conversion function to Salinity, Resistivity and TDS



Dissolved Oxygen Meter OM-71

DO	Saturated Oxygen	Oxygen
----	------------------	--------

- Compact design with special shock and alcohol resistant body
- Fully waterproof/dustproof meter and electrode (IP67 rated)
- Large LCD with backlight
- Data memory up to 1000 sets and PC/Printer output
- Air calibration, Salinity concentration correction



Multiparameter Water Quality Checker U-50 series

pH	ORP	COND	SAL	TDS	DO	TURB
DEPTH	TEMP	GPS				


- Ideal for water quality testing and inspection of river, lake, well water, groundwater, discharge water and other water sources
- Simultaneous measurement and display of up to 11 parameters
- Integrated sensor probe and display section for maximum portability
Convenient for one-point measurement and measurements near the surface of the water. Built-in highly sensitive turbidity sensor enables measurement of even low turbidity water



Water Quality Monitoring System W-20XD Series

pH	ORP	COND	SAL	TDS	DO	TURB
DEPTH	TEMP	ION				

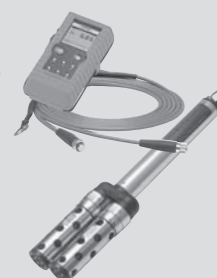- Ideal for water quality testing and inspection of, city sewage water, lakes and marshes, dams, wells and ground water, factory drainage, farm water, and nurseries
- Simultaneous measurement and display of up to 13 parameters
- Up to one month data logging (With measurements every 15 minutes)
- Measurement at depths as low as 100 meters with its superior durability and high pressure resistance



Compact LAQUAtwin

pH METER B-711/B-712/B-713

pH

- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.1 mL
- Select measurement method depending on your situation and sample. (Drops, Immersion, Scoop, Wipe, Solid samples, Powders, Paper, textiles)
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold

B-711 (One-point calibration)

B-712 (Two-point calibration)

B-713 (US only) (Two-point calibration)



pH

Compact Conductivity METER B-771

COND

Salt

TDS

- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.12 mL
- Conductivity readings can be converted into Salt and TDS
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change



COND

Sodium Ion METER B-722

ION

- Only compact meter for a quick measurement of sodium ion using ion selective membrane
- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.3 mL
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change



Na⁺

Potassium Ion METER B-731

ION

- Only compact meter for a quick measurement of potassium ion using ion selective membrane
- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.3 mL
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change



K⁺

Nitrate Ion METER B-741/B-742/B-743

ION

- Only compact meter for a quick measurement of Nitrate ion using ion selective membrane
- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.3 mL
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change
- Special application packages for crop and soil

B-741 (for crops) B-742 (for soil)

B-743 (for general use)



NO₃⁻

Calcium Ion METER B-751

ION

- Only compact meter for a quick measurement of Calcium ion using ion selective membrane
- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.3 mL
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change



Ca²⁺

Salt METER B-751

ION

- Only compact meter to measure sodium ion to calculate into NaCl based salt concentration unlike the conductivity converted meters
- Flat sensor technology realizes a reliable and direct measurement of a drop of the sample from 0.3 mL
- IP67 waterproof and dustproof
- Temperature compensation/Auto hold/Auto range change



Salt

WATER QUALITY ANALYZERS website

With more than 60 years of engineering excellence, HORIBA's diverse lineup of water quality analyzers and electrodes can be easily suited for any laboratory's usage. Our Water Quality Analyzer website provides you with a wide range of information that you might find useful; such as water quality measurement tips and downloading of manuals for our meters and electrodes.

<http://www.horiba-water.com>



Benchtop Meter (LAQUA) website

Taking your comments and feedback to our hearts, "LAQUA" is our new brand which delivers you the best solution for your needs. Our "LAQUA" website features a "Selection Guide" which enables you to find, within HORIBA's variety of meters and innovative electrodes, the most suitable one for your usage.

<http://www.horiba.com/laqua>



Compact Meter (LAQUAtwin) website

Analyzing water quality will be more simple for you when using our LAQUAtwin range of meters. Anyone, anywhere, at any time, You will be able to measure samples easily with the LAQUAtwin. You can experience the charm of the LAQUAtwin at our LAQUAtwin website.

<http://www.horiba.com/laquatwin>



Please read the operation manual before using this product to assure safe and proper handling of the product.

- The contents of this catalog are subject to change without prior notice, and without any subsequent liability to this company.
- The color of the actual products may differ from the color pictured in this catalog due to printing limitations.
- It is strictly forbidden to copy the content of this catalog in part or in full.
- All brand names, product names and service names in this catalog are trademarks or registered trademarks of their respective companies.